



CC-CLIP


# CLIP FEATURES

**1.** The Calvicular Clip loop hides the interconnecting suture knot into the clavicular drill hole: eliminates protruding knots.

**2.** The outer diameter 5.0 mm, is a snug fit into the 5.5-6.0-mm clavicular drill hole. It squeezes the tendon graft against the drill hole walls: eliminates the use of an interference screw.

**3.** The width (8.0 mm) and length (16.0-20.0 mm) of the Calvicular Clip provides a snowshoe effect: eliminates implant sinking into the bone.

**4.** The dorsal rim of the Calvicular Clip reduces pressure induced by the tendon graft sling onto the clavicular bone: eliminates bone resorption.

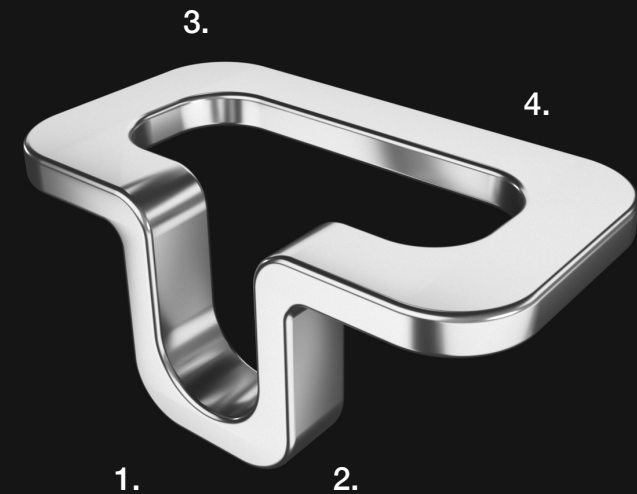


# CC-CLIP PRODUCTS

REF. NUMBER

JR-002

JR-003

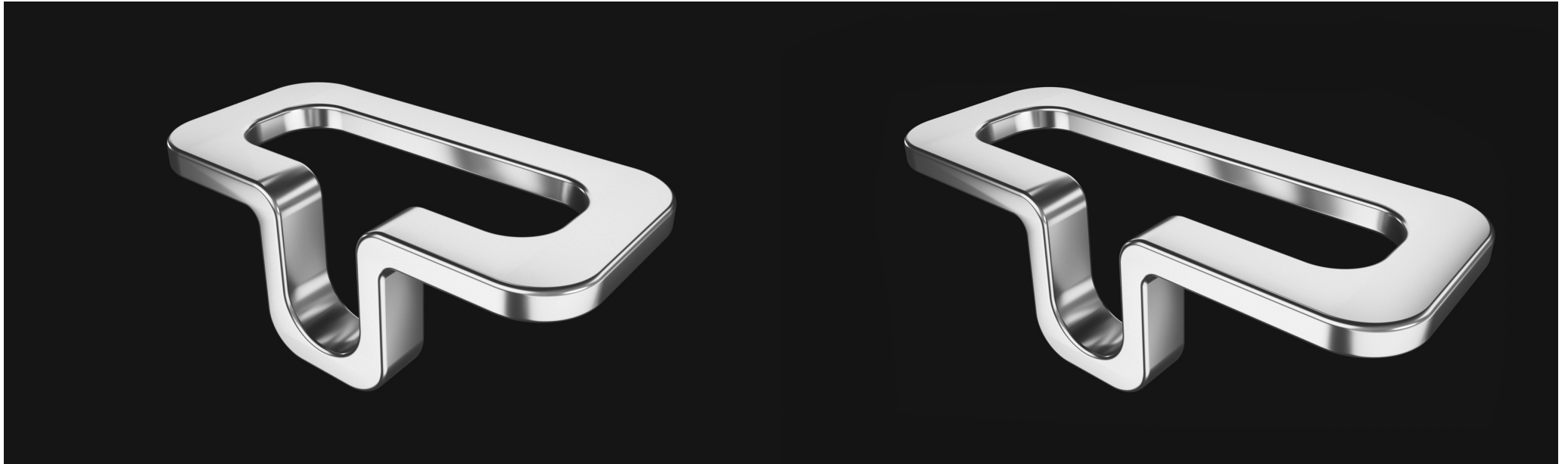

DESCRIPTION

Clavicle Clip

16 mm

Clavicle Clip

20 mm

JR-007

JR-008

JR-011

JR-012



Sub-Coracoid Clip  
10 mm


Sub-Coracoid Clip  
12 mm

Sub-Coracoid Clip  
10 mm, round

Sub-Coracoid Clip  
12 mm, round

JR-023

JR-020

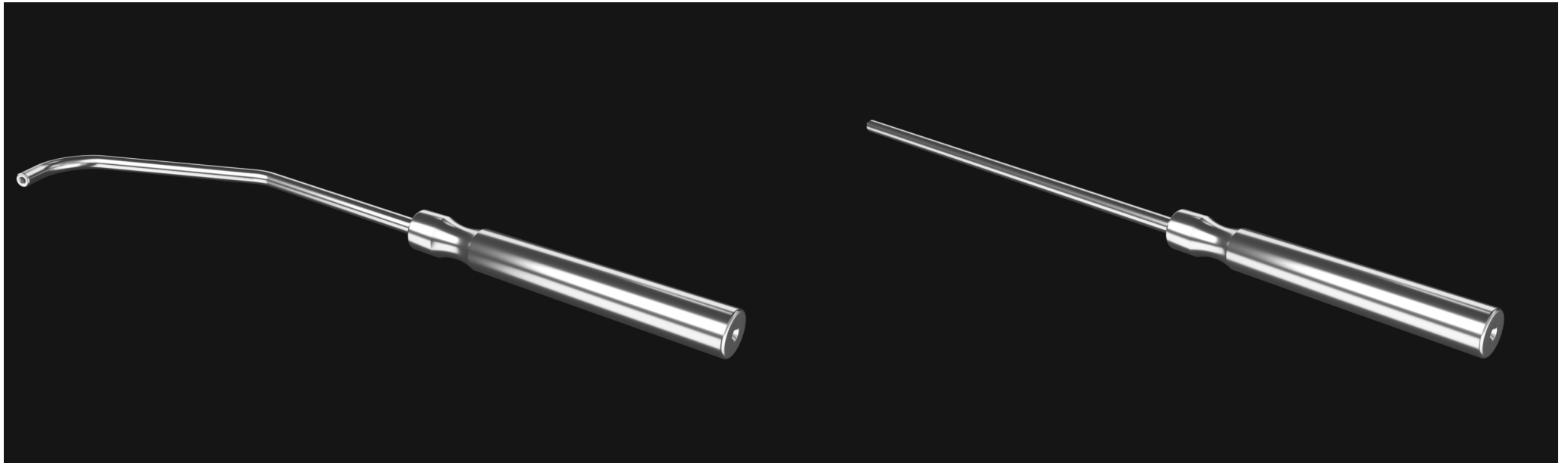

Lasso Guide Curved

D7,30xd1,10 L=165

Lasso Guide Straight

D7,30xd1,10 L=165

